


Aesculap® Safety scalpels

Certified acc. TRBA250

CE0123

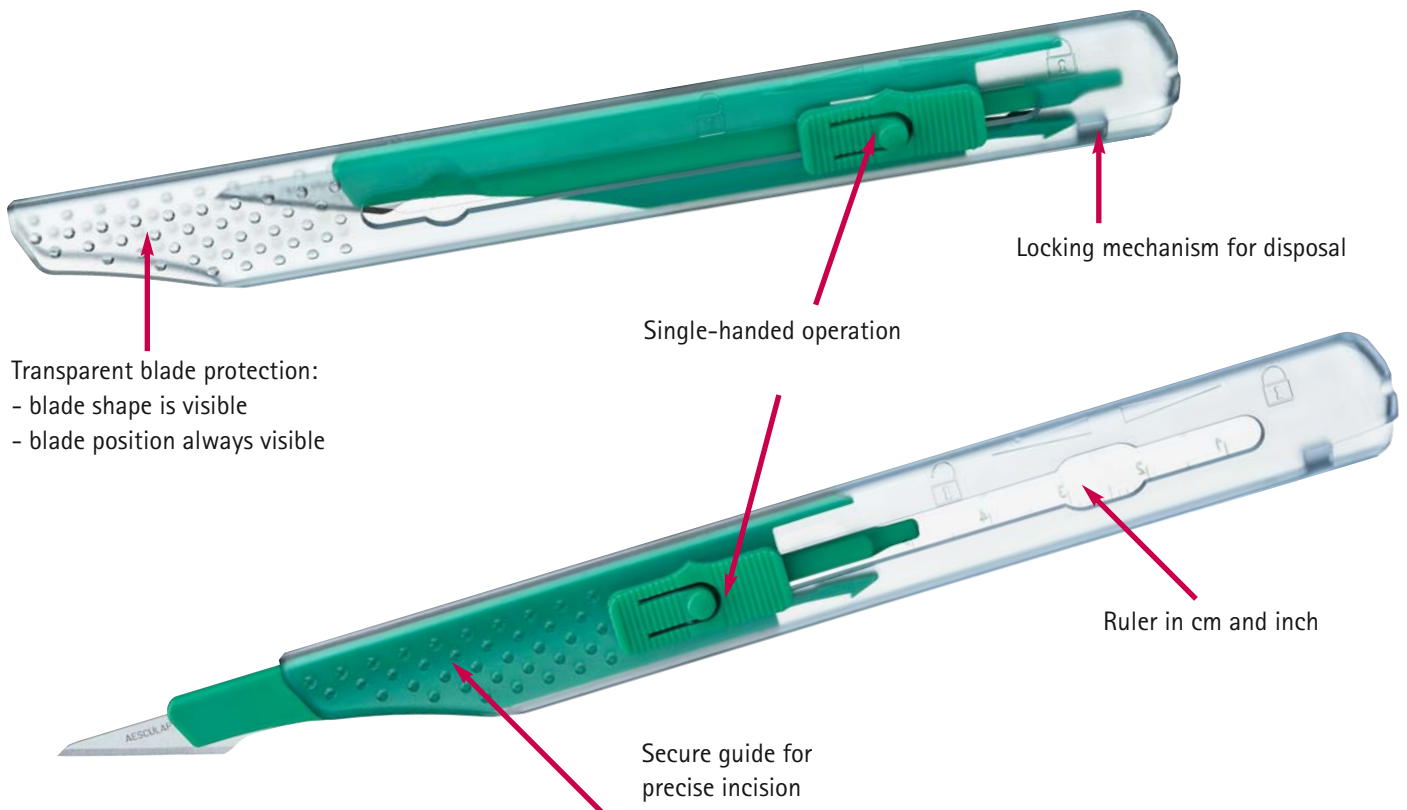

Transparent blade protection:
- blade shape is visible
- blade position always visible

Single-handed operation

Locking mechanism for disposal

Ruler in cm and inch


Secure guide for precise incision



For further information see our brochure C65402.



Comparison between the new TRBA 250 requirements and features of the Aesculap safety scalpel

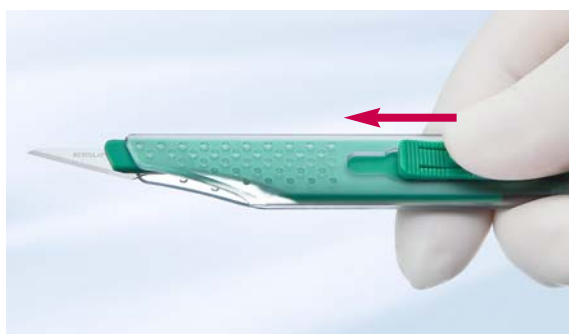


§ The safety mechanism is part of the system and compatible with other accessories

§ Single-handed activation must be possible

✓ The mechanism is activated automatically by sliding back the blade with one hand

§ The safety mechanism must prevent reuse

✓ The blade is safely locked in the disposal position

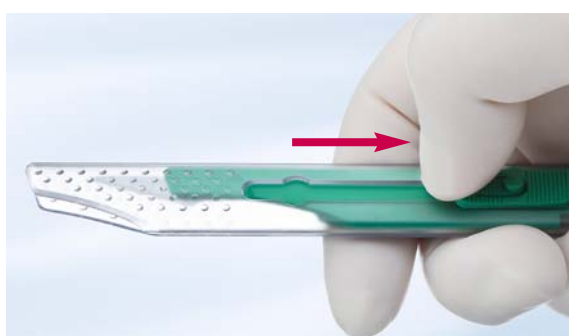

§ Using the safety product does not require a different technique of use.

✓ It is used in the same way as a standard single-use scalpel

§ The activation of the safety mechanism must be clearly indicated (palpable or audible)

✓ An audible "click" indicates that the safety locking mechanism has been activated, and this can also be seen immediately



"click"

palpable

Aesculap® Safety scalpels

Certified acc. TRBA250

CE0123



BA810SU-BA836SU Sales unit PAK = Package of 10 pieces



BA819SU
Fig. 19



BA820SU
Fig. 20



BA821SU
Fig. 21



BA822SU
Fig. 22



BA823SU
Fig. 23



BA824SU
Fig. 24



BA825SU
Fig. 25



BA836SU
Fig. 36