

# AquaCare

Contactless Dentistry )))

*Cleans, Cuts, Prepares - Gently*

## Contactless Dentistry •

- Unlike conventional rotary cutting instruments, the AquaCare handpiece is not in direct contact with the tooth structure, removing only the minimum of sound tooth material and **eliminating the risk of chipping and stress fractures**.
- AquaCare **increases bond strength** with its texturing of the cavity surface.
- **Faster, more efficient cleaning** for removing stains from all tooth aspects, even accessibly difficult pits and crevices, whilst also treating sensitivity. Patients love their brilliant smiles.



- Minimal need for anaesthetics combined with the lack of vibration, heat generation and noise leads to **reduced patient anxiety**. AquaCare is recognised as a practice enhancing product.
- Balanced, slim line **fully autoclavable handpiece** without bulky attachments.
- Incredibly **versatile** across a wide variety of restorative dental procedures.
- Convenient, **easy to use** colour-coded cartridge system makes swapping between types of treatment effortless.




# AquaCare Twin

## Superior polishing and prophylaxis

ProSylc powered by NovaMin®

- Treatment of dental hypersensitivity.
- Fast removal of extrinsic stains, e.g. coffee, tea, tobacco.
- Plaque removal prior to fluoride treatments.
- Pre and post whitening treatments.

CLEANING

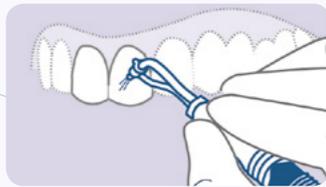

Angle 45°. Approx. 4mm  
away from the surface.

Use AquaCare before scaling to remove stubborn plaque and decay quickly, leaving only very tough calculus for removal with a scaler.

## ProClean for Sodium Bicarbonate

- Fast stain removal with no patient discomfort.




## Superior Abrasion & Preparation

ProCut 29µ & ProCut+ 53µ Aluminium Oxide

- Roughening of surface in preparation for adhesive dentistry.
- Prevention of micro-fractures and crazing normally associated with rotary instrumentation.
- Opening and clearing of fissures.
- Enhanced sealant retention.
- Removal or repair of composite restorations.
- Cavity preparations for composite restorations and small amalgams.
- Preparation of small carious lesions.

BONDING / PREPARATION

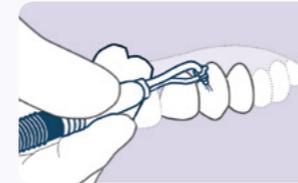

Angle 35-45°. Approx. 2-3 mm  
away from the surface.

BONDING / PREPARATION




Angle 45°. Approx. 2mm  
away from the surface.

- Matting the surface of metallic restorations for detection of interferences.
- Remove cement from crowns, bridges and brackets.

## It's a Treatment Centre

- Designed to deliver abrasion and prophylaxis, the AquaCare Twin Treatment Centre empowers the clinician to switch easily between dental procedures.
- Pressure and media settings do not need adjusting for differing treatments.
- Enjoy the finesse required in multi-disciplinary treatments.

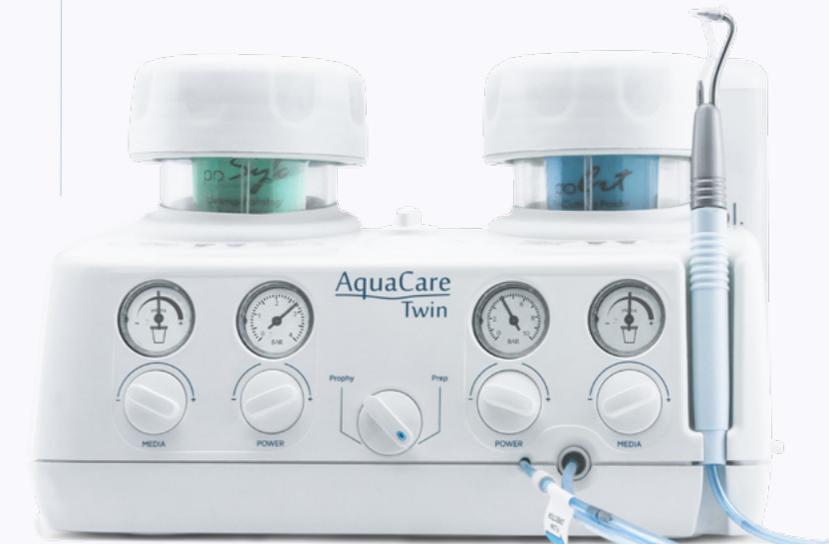




**4-in-1 Handpiece**  
cut | clean | wash | dry



- 'Venturi' design blockage free tip.
- Intuitive, easy-to-use power, fluid and media controls.
- Ultra-lightweight [less than 30g] fully autoclavable handpiece.
- 0.6mm and 0.8mm handpiece options.
- Single patient disposable nozzle and tube.
- Intelligent 3-position foot control delivering:

air | air + solution | air + solution + power



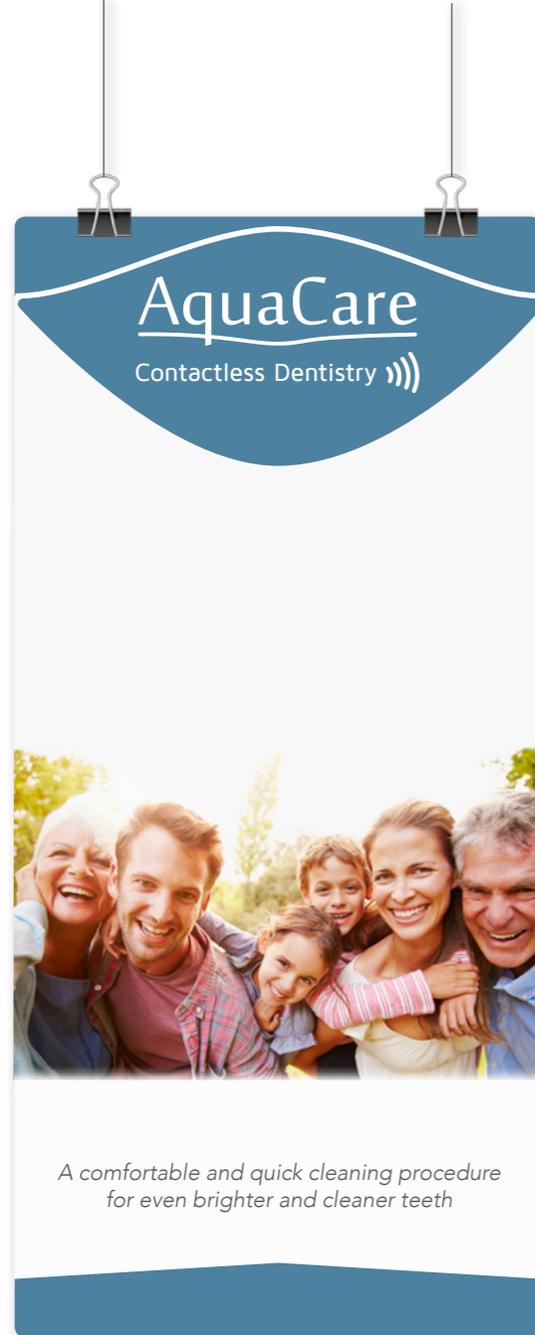

## AquaCare

AquaCare uses breakthrough technology to deliver a comfortable, quick cleaning procedure by using a fine stream of fluid combined with a tiny volume of powder directed at the teeth to be treated. Decay is then comfortably removed with minimal sensation, leaving healthy tooth material undisturbed.

**AquaCare Patient Leaflet** designed to make it easy for you to explain the benefits of air polishing to your patients. Find out more at [velopex.com](http://velopex.com) or call us on +44 (0)20 8965 2913



[www.velopex.com](http://www.velopex.com)



Medivance Instruments Ltd., Barretts Green Road,  
London NW10 7AP United Kingdom  
+44 (0)20 8965 2913 [enquiries@velopex.com](mailto:enquiries@velopex.com)