



# FUSION<sup>5</sup> Curing Light

*The Most Powerful and Dependable Curing Light Platform!*

*3 seconds cure and much more ...*

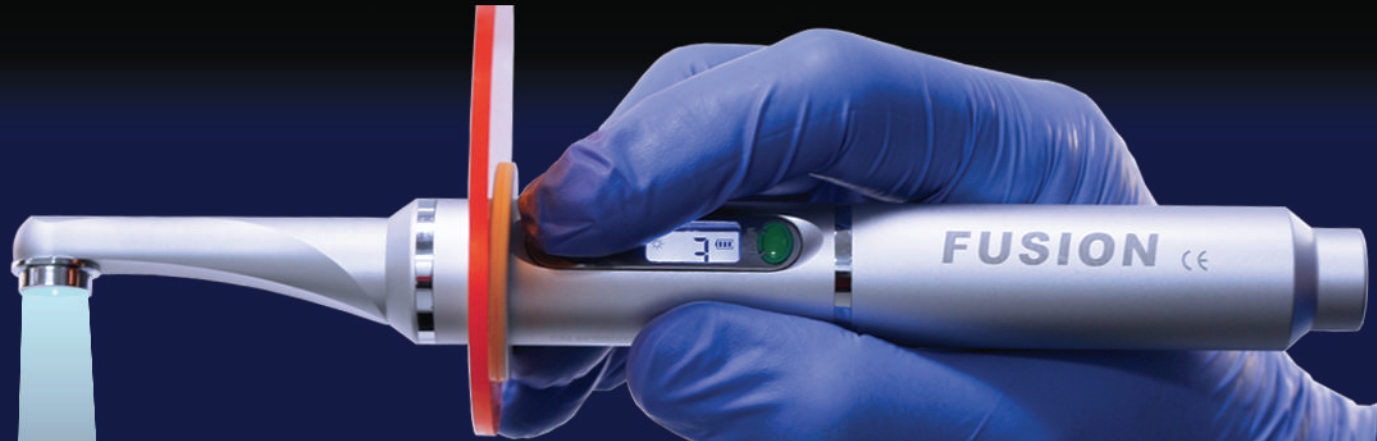

- 3-second cure with focused beam
- Flexible platform for oral exams

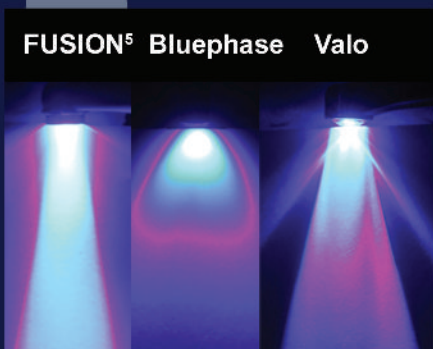

FUSION<sup>5</sup> Bluephase Valo

Focused      Dispersed



FUSION Curing Light is a 21st century light that is versatile, cost-effective, and valuable in everyday practice. Adding to the value is the ability to conduct oral cancer examination.

**Dr. Howard Glazer**



The Fusion Light is one GREAT all around light for the dental office. You get the benefits of three lights in one! Curing, Transillumination and Oral Cancer Screening. The Curing power is remarkable.

**Dr. John Comisi**



Version 5

**DentLight**

<https://stomshop.pro>




# FUSION<sup>5</sup>

## Amazing Power and Style in the Palm of Your Hand!

- Unique focused beam and advanced backlit LCD display:
  - 4000\* mW/cm<sup>2</sup> Plasma Power mode for fast/deep cure
  - 2000\* mW/cm<sup>2</sup> Pulsed mode for low heat/shrinkage cure
  - Battery and thermal bars for smart control
- Wide spectrum to cure all materials
  - 420 - 490nm
- Sleek metal keypad with large activation and program buttons
- Modular design with easy upgrade and service parts
- Cordless with direct power plug that eliminates downtime
- 90 minutes fast recharge time
- Super slim and ergonomic
  - Tip height only 11.5mm
  - Weighs only 3.9oz
  - Measures at 7.8 inch

(\* Absolute Irradiance using Optical Spectrometer with Cosine Corrector Detector)

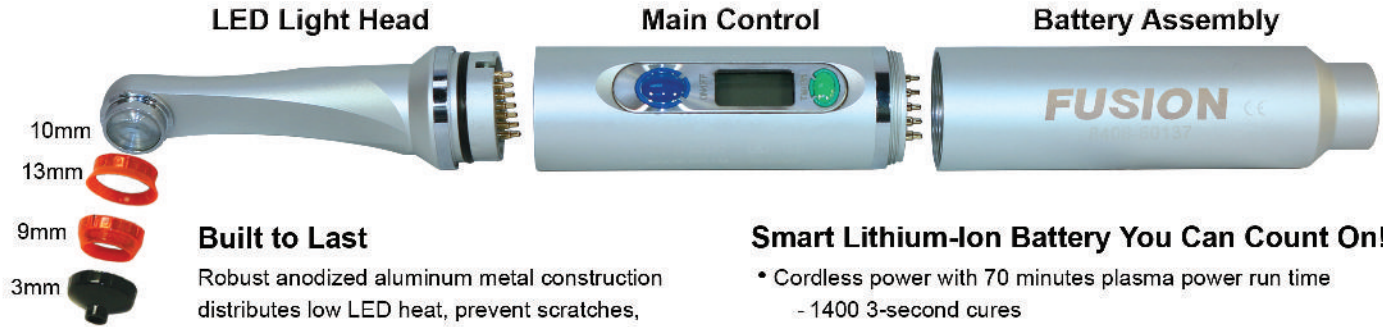

### Built to Last



Robust anodized aluminum metal construction distributes low LED heat, prevent scratches, discoloration and damage from accidental falls.

### Smart Lithium-Ion Battery You Can Count On!


- Cordless power with 70 minutes plasma power run time
  - 1400 3-second cures
- Constant output throughout battery run time
- Common platform for future upgrades

### Lowest Profile

Easily cure posterior teeth with small 90° head

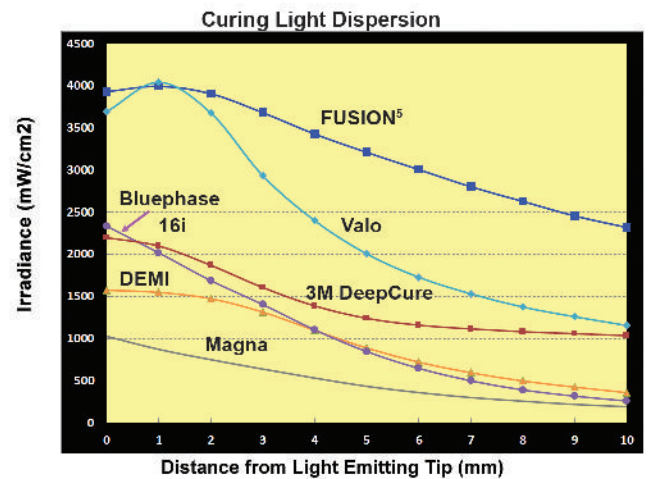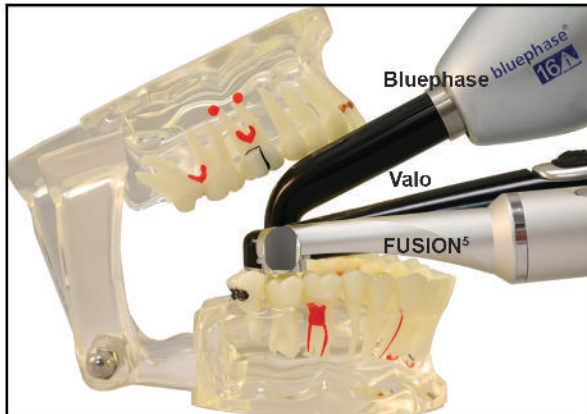

### Oral Exam Functions with Interchangeable Heads

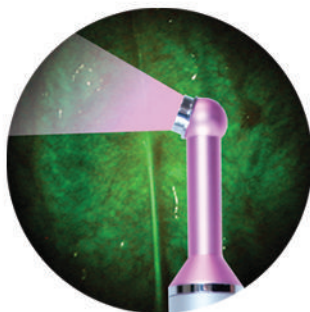

Oral Exam (Cancer/Caries Detection)  
PN 7500152



Cracks/Canal Detection  
7800060W

# PROFIT FROM the Award-winning FUSION platform

## Order Information

### PN:7800080 FUSION<sup>®</sup> Curing Light Kit:

Including Standard Handpiece, Charging Stand, Curing Caps (tacking, 9mm and 13mm), Power Adapter, Light Shield, Barrier Sleeves (100), and Instructions for Use.



### Combination Curing and Oral Exam Options

#### 7700512 FUSION-DOE Kit

Including FUSION<sup>®</sup>, Violet Light Head, Green and Orange Fluorescence Goggles, Curing Caps, Charging Stand, Power Adapter, Barrier Sleeves, and DVD or USB training guide, and Metal Carrying Box.

#### 7710010 FUSION-DOE Deluxe Kit

Including FUSION<sup>®</sup>, Violet and White Light Heads, Green and Orange Fluorescence Goggles, Curing Caps, Charging Stand, Power Adapter, Barrier Sleeves,(100), and DVD or USB training guide, and Metal Carrying Box.

### Add-on Options and Accessories

7500152 FUSION<sup>®</sup>-DOE starter Kit

8500123 Green Fluorescence Goggle

7500088 Orange Fluorescence Goggle

7700192 Orange Fluorescence Shield



7700092 Green Fluorescence Shield

7500060U Violet Light Head


7500082 Barrier Sleeves (100/pack)

7701100 Green Camera Filter (fit on SLR)

7701110 Orange Camera Filter (fit on SLR)

8500122 Green Fluorescence Loupe Filter Pairs

8501122 Orange Fluorescence Loupe Filter Pairs

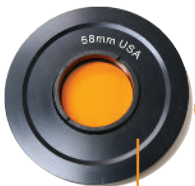

# DentLight

<https://stomshop.pro>

1825 Summit Ave. #210

Plano TX 75074

Tel: 972-889-8857

Email: [sales@dentlight.com](mailto:sales@dentlight.com)

[www.dentlight.com](http://www.dentlight.com)