

E-FLOW

Dental anesthesia delivery system

End / Pain / Now



Patented Hydraulic Drive Technology

Painless Performance




*High-precision sensor is used for more precise control, the injection pressure is just hydraulic pressure.

<https://stomshop.pro/>




4.3 " Touch Screen

Visualization of various functions, bringing sensory experience of visual interaction

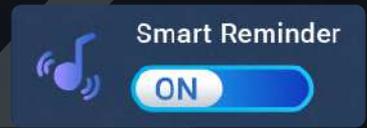

Smart Voice Reminder

Timely voice reminders under low anesthesia dose and high system pressure(in PDL mode), which is safer and more secure

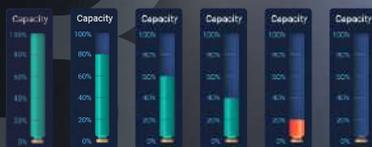

High-precision Liquid Control

Not only can achieve precise drug delivery accuracy, but also can see the real-time remaining dose



Real-time Dynamic Pressure

The real-time injection pressure is displayed through the pressure value and pressure histogram, bringing a painless and safe treatment

HANDPIECE



Autoclavable Container

Prevent cross-infection, safe and economical



Aspiration

Prevent swelling and pain caused by excessive liquid

Lighter and Smaller

Only weighs 40g, easier to operate

Intelligent Foot Pedal Control

The injection flow rate can be switched by pressing the foot pedal to the "front section" or "back section" of the stroke



FOOT PEDAL

4 injection modes, multiple options for effective pain relief



PDL Mode (periodontal ligament)

Automatic injection can be selected in this mode

Dosage Mode

Free choice of injection dose,
making the injection more stable and precise





| | |
|-------------------|--|
| Model: | E-FLOW |
| Weight: | 1.8kg±10% |
| Handpiece weight: | 40g |
| Injection mode: | Normal, Fast, PDL, Dosage |
| Injection speed: | Slow (0.3ml/min) Medium (1.7ml/min) Fast (3.5ml/min) |
| Battery capacity: | 2600mAh |

SPECIFICATIONS